




Linear Cutter Stapler

SSAB




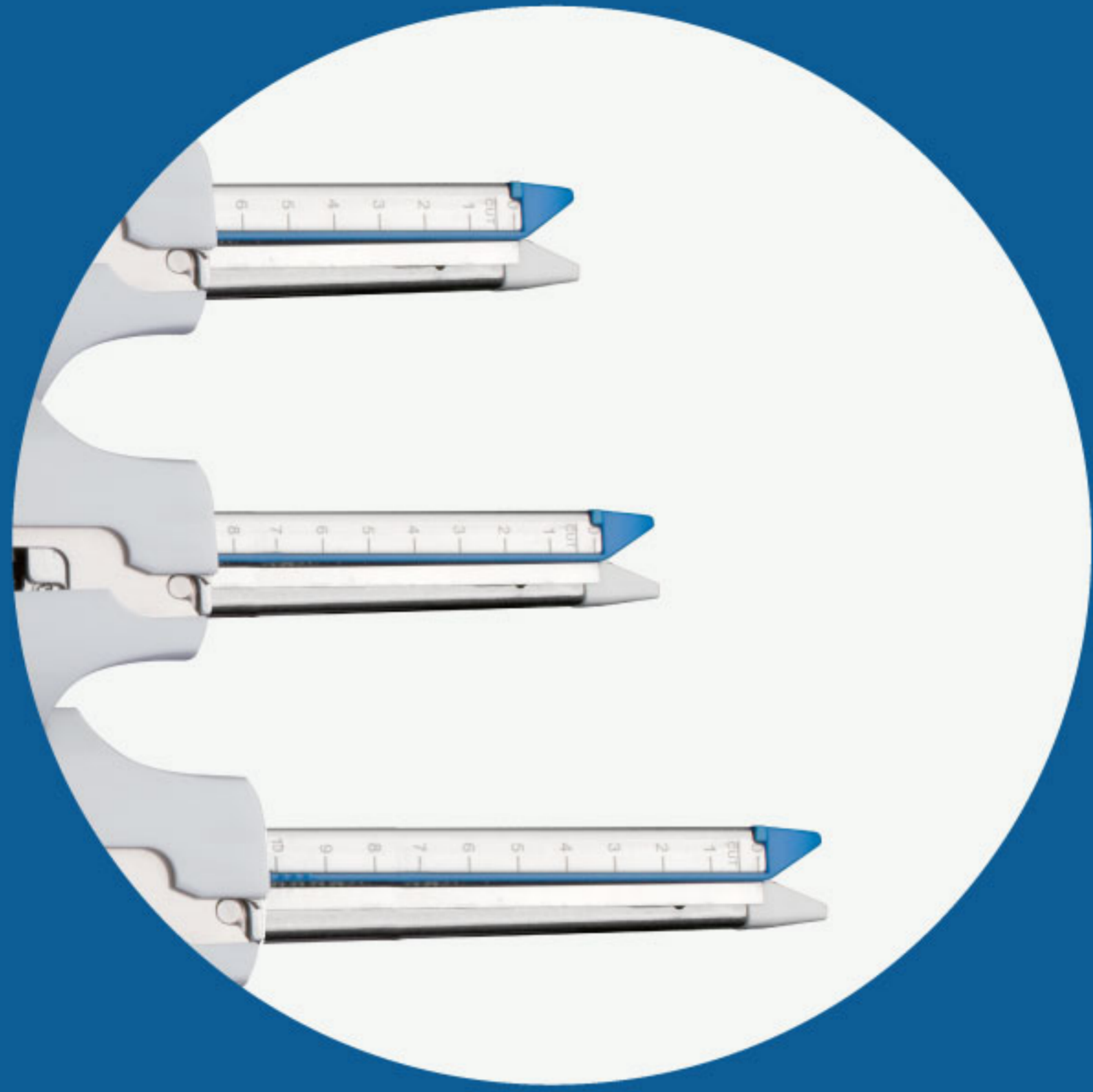


Three different staple line lengths of reloading unit, 60 mm, 80 mm and 100 mm, addressed to various tissues



Locking handle design, can be used in one hand operation

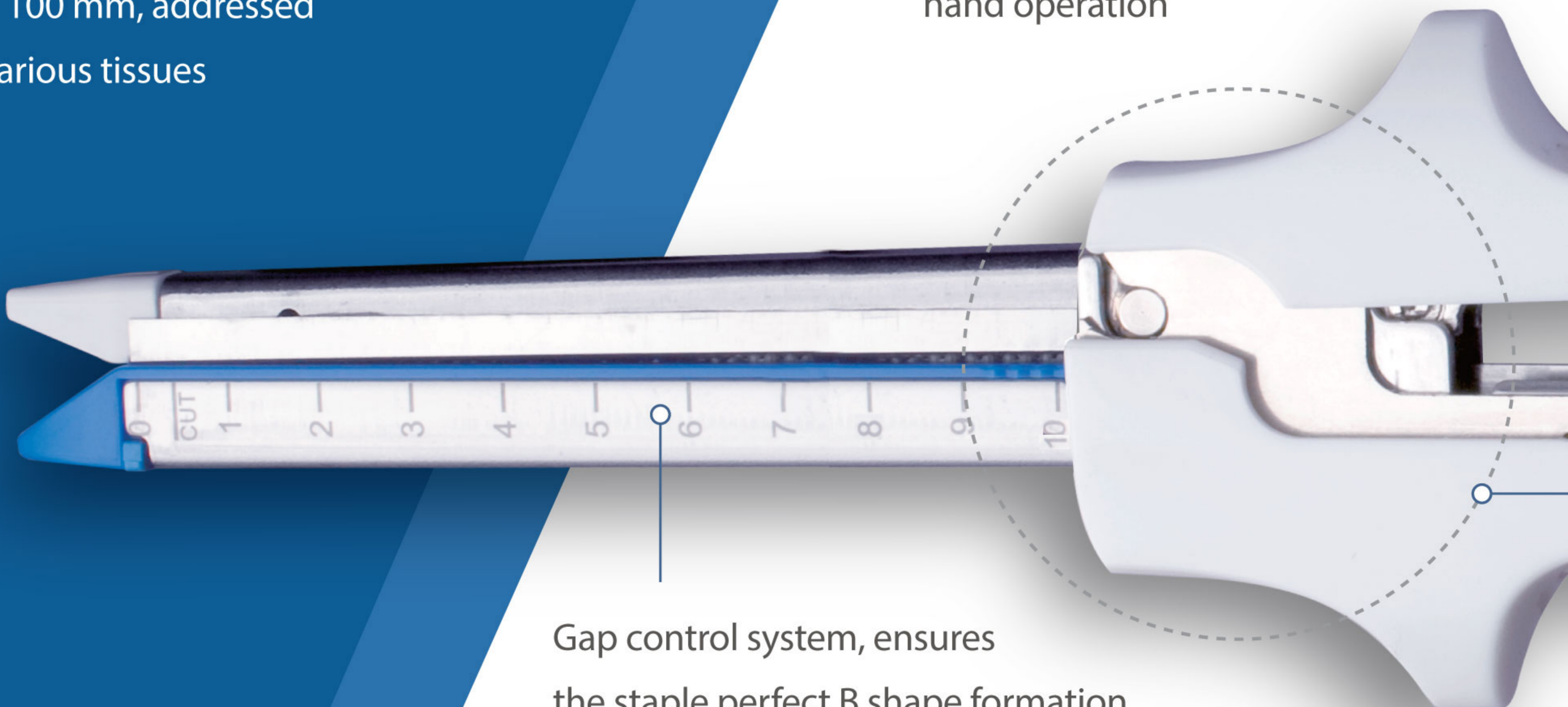


Gap control system, ensures the staple perfect B shape formation from proximal to the distal end

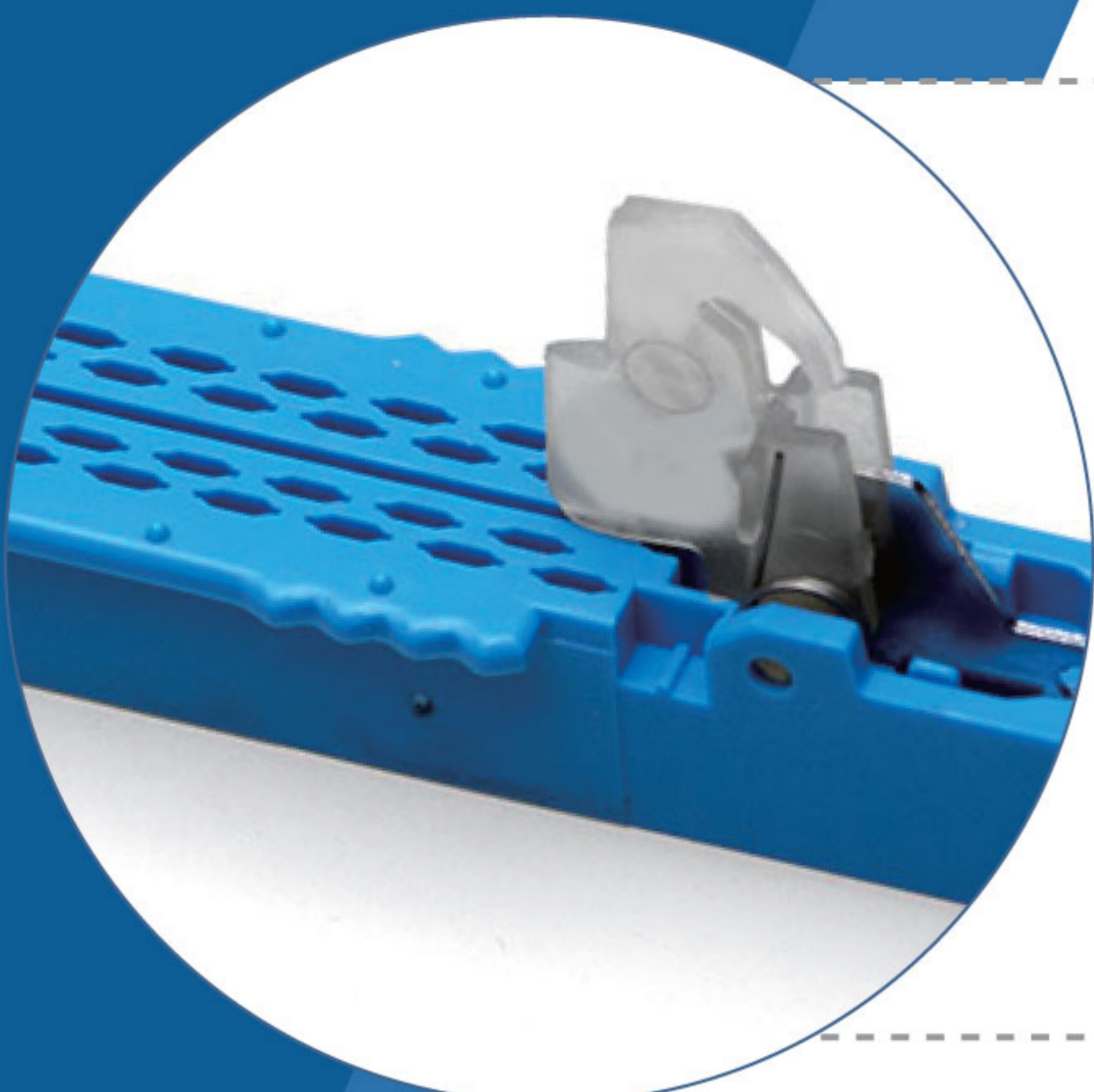

Blade protector and cartridge safety lockout:
White plastic protector covers blade before and after reloading unit fired, which prevents surgeon or assistant cutting finger on reloading unit blade.
Authentic disposable blade, ensures effective cutting, reduces damage on tissue, avoids cross contamination


SSAB LINEAR CUTTER STAPLER

Panther SSAB

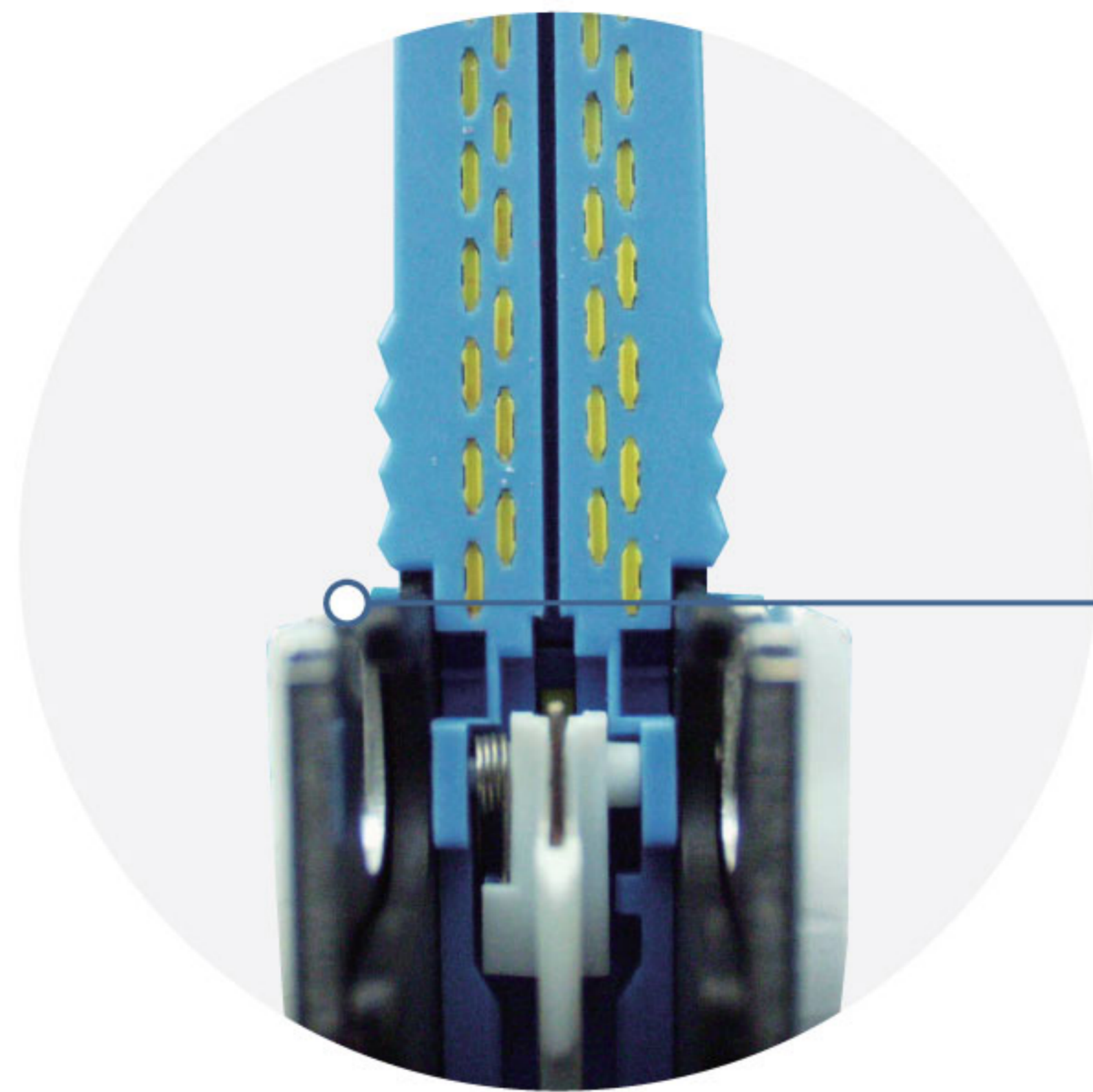


Tissue end position

Other brands

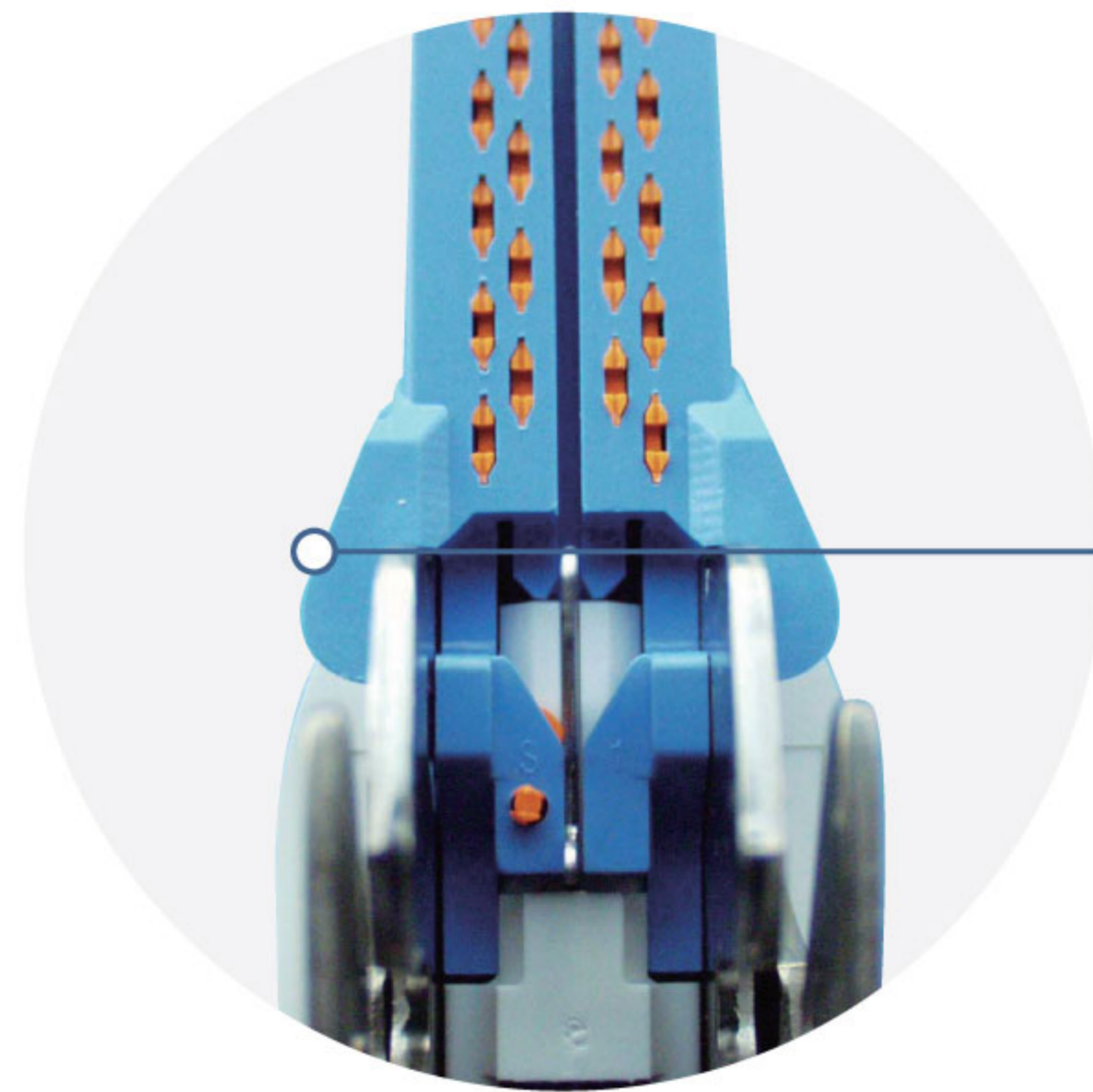


Tissue end position

Tissue end position, ensures the tissue to be effectively stapled after cutting, much safer

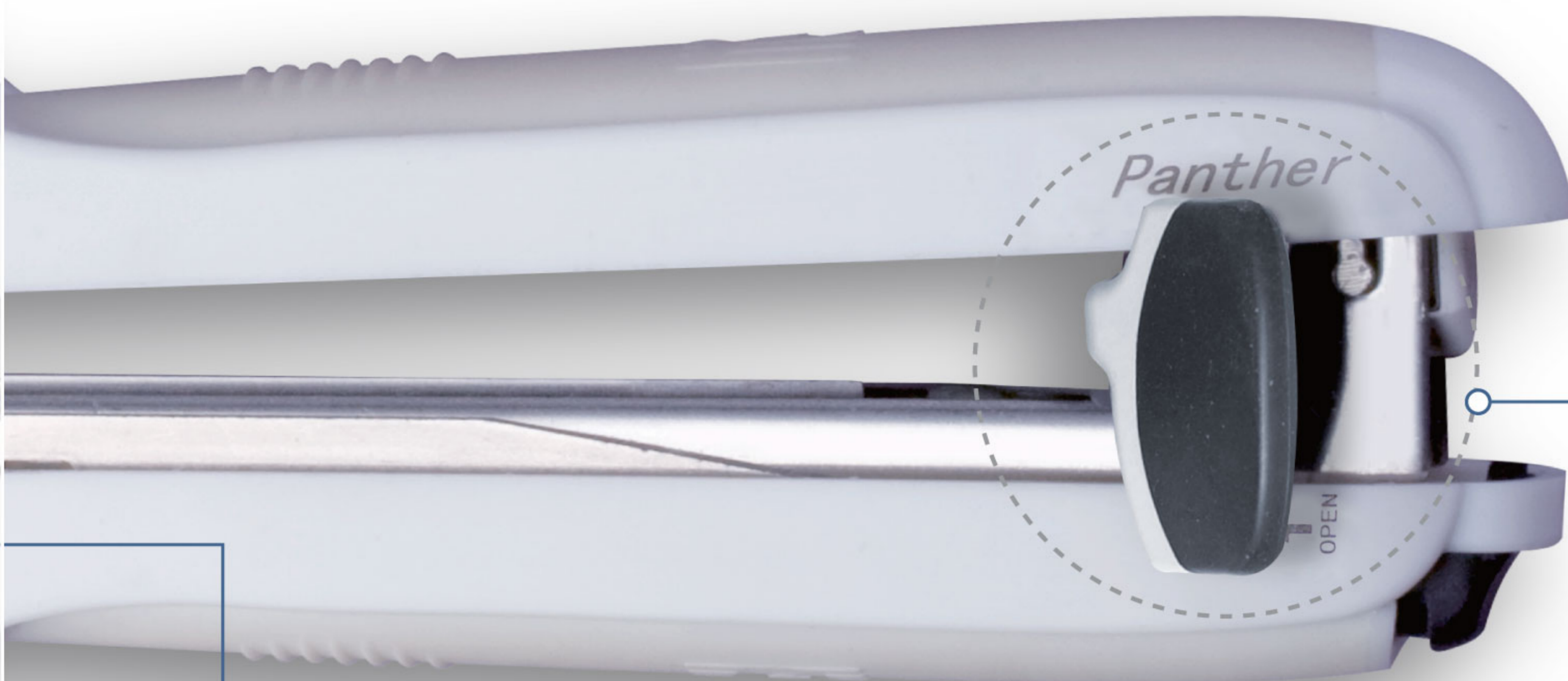

Adjustable firing knob:
Ability to fire from either side of stapler, flexible and easy to use

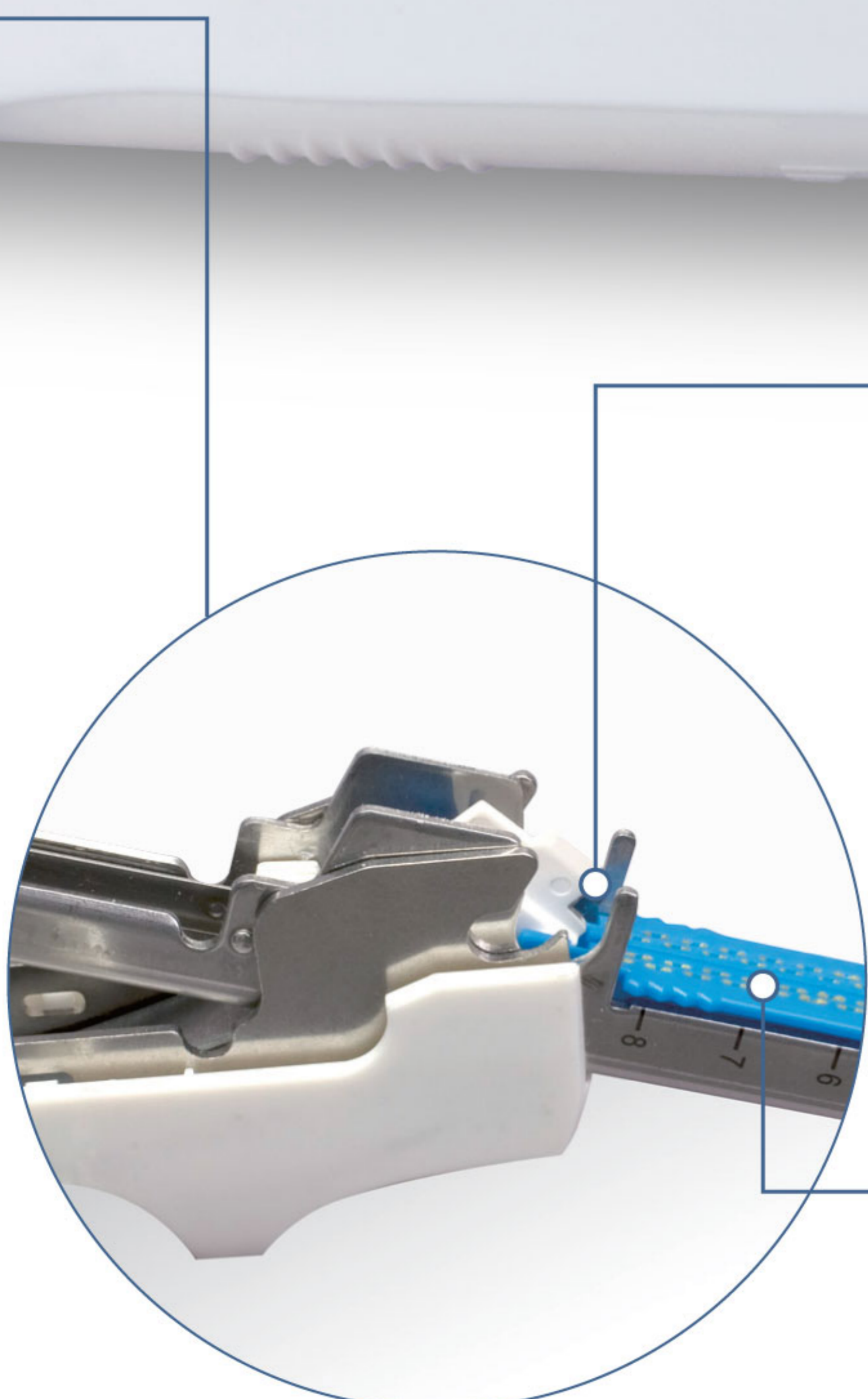

New blade with each cartridge:
Guarantees effectively cutting ,
reduces tissue trauma, avoids
cross contamination

Floating reloading unit design:
The reloading unit floats above
the edge of the fork before
closing the stapler fully,
avoids mis-firing

LINEAR CUTTER STAPLER AND RELOADING UNITS



Product Name	Product	Description
Linear Cutter Stapler	SSAB-60 (3.8mm)	Linear Cutter Stapler 60mm, 3.8mm staple height
	SSAB-60 (4.8mm)	Linear Cutter Stapler 60mm, 4.8mm staple height
	SSAB-80 (3.8mm)	Linear Cutter Stapler 80mm, 3.8mm staple height
	SSAB-80 (4.8mm)	Linear Cutter Stapler 80mm, 4.8mm staple height
	SSAB-100 (3.8mm)	Linear Cutter Stapler 100mm, 3.8mm staple height
	SSAB-100 (4.8mm)	Linear Cutter Stapler 100mm, 4.8mm staple height
Linear Cutter Stapler Reloading Unit	SADB-60N	Linear Cutter Stapler Reloading Unit 60 mm, 3.8mm staple height
	SADB-60D	Linear Cutter Stapler Reloading Unit 60 mm, 4.8mm staple height
	SADB-80N	Linear Cutter Stapler Reloading Unit 80 mm, 3.8mm staple height
	SADB-80D	Linear Cutter Stapler Reloading Unit 80 mm, 4.8mm staple height
	SADB-100N	Linear Cutter Stapler Reloading Unit 100 mm, 3.8mm staple height
	SADB-100D	Linear Cutter Stapler Reloading Unit 100 mm, 4.8mm staple height